



Sampling Instructions - DOG

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

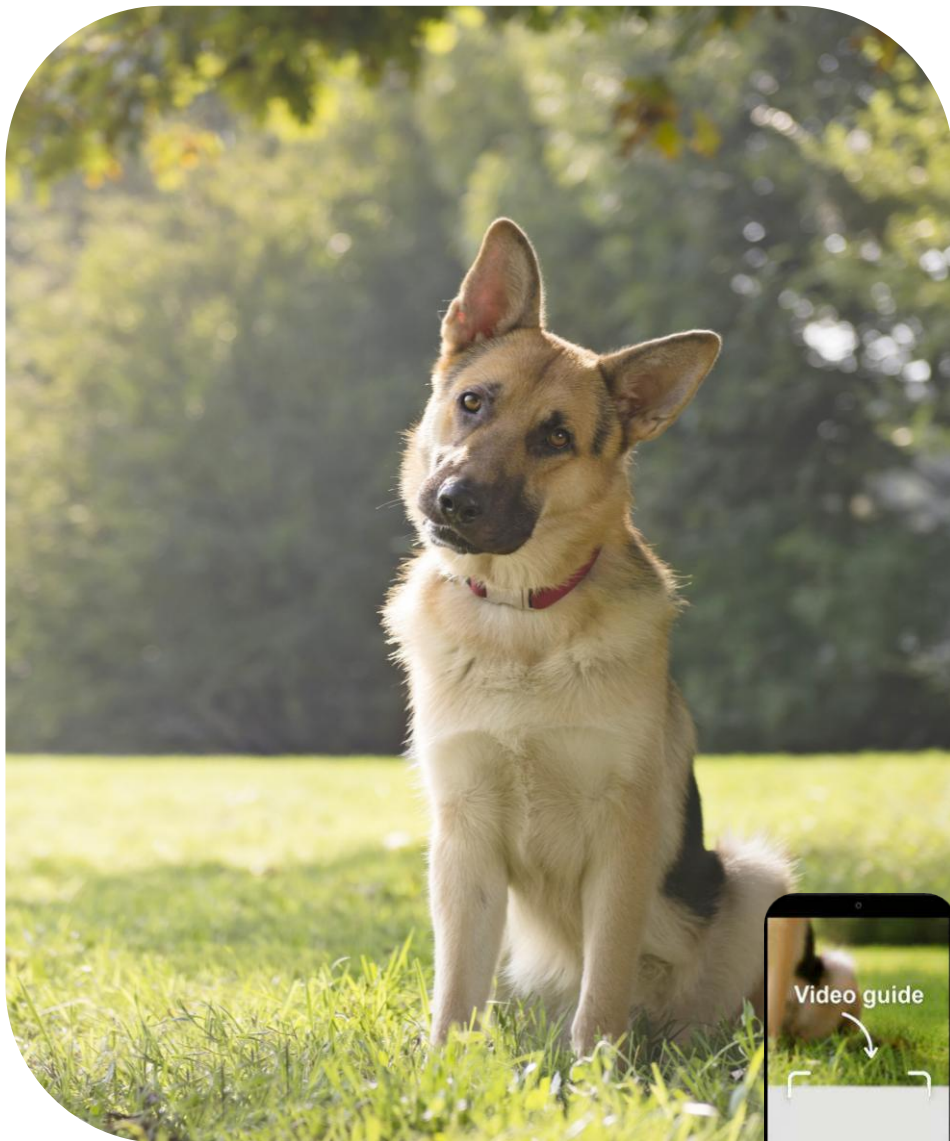

- ✓ **GOcard/FTAcad blood card**
- ✓ Mucous membrane swab
- ✓ EDTA blood sample (vet. staff only)

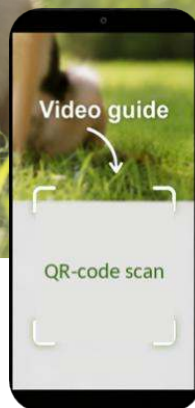

Sampling Instructions - DOG

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com


GOcard/FTAcad blood card

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organization, follow the applicable regulations


Preparation

□ Fill in GOcard / FTAcad

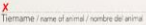GOcard®-SC Tierregistrierung Animal registration Inscripción del animal generatio Center for Animal Genetics & Animal Trust Center	Generatio GmbH Center for Animal Genetics & Animal Trust Center Blumenstr. 49 69115 Heidelberg Germany - Alemania	 Die Nummer des Online-Auftrags lautet: Associated online order no. La correspondiente orden en línea # X Ein Auftragsformular liegt bei: Order form is attached Se adjunta una lista de los exámenes
 X Einsender / sender / remitente	 X Tiername / name of animal / nombre del animal	 X Chip / tattoo
 X Unterschrift / signature / firma		


Sampling

Clean ear (alcohol swab)

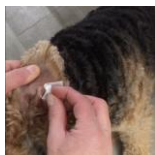

1

Remove lancet safety pin by twisting

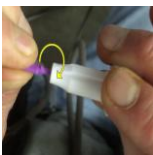

2

Place lancet firmly, then release



3

Wait for blood drop, massage towards puncture site if needed




4

Apply blood to clean GOcard



5

Allow blood drop to air-dry

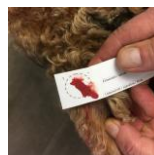

6

Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax

Sampling Instructions - DOG

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

GOcard/FTAcad blood card


Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organization, follow the applicable regulations

Loaded GOcard/FTAcad

	X Einsender / sender / remitente	X Tiername / name of animal / nombre del animal	Notes
	X Unterschrift / signature / firma	Chip / tattoo	

Prepare sample shipment

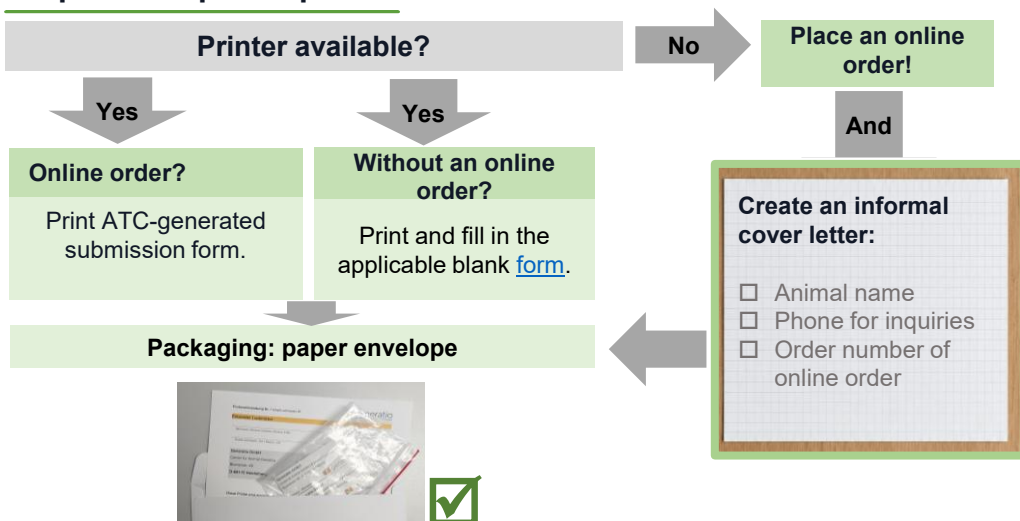

Shipping instructions



Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax

Sampling Instructions - DOG


Mucous membrane swab

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organization, follow the applicable regulations

Preparation

☐ Label tube for mucous membrane swab

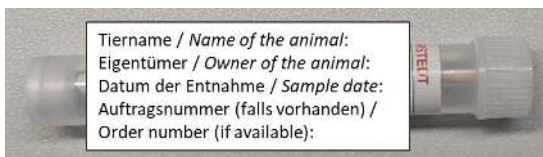

Sampling

- 1 Open tube by twisting and pull out the swab stick.

 - 2 Carefully insert the swab between cheek and gum.

 - 3 Rub the swab firmly against the inside of the cheek for approx. 20 seconds.

 - 4 Place swab directly back into tube

- Repeat with 2nd swab

 With our collection kits, the swabs don't have to be dry when you put them back in the tube, saving you 40 minutes.

Notes

-  at least 2 hours after the last meal.
-  Do not let the swab contact anything other than the sample.
-  Separate animal from others — exclude DNA contamination.
-  For puppies: take swab only after weaning.

 If you are sending us mucosal samples without a collection kit, please use standard cotton swabs. Please allow the swabs to air dry for at least 20 minutes. Please ensure that swabs from different animals do not come into contact with one another (for example, by placing the swabs individually in paper envelopes and labeling each one according to the sample preparation instructions).

Sampling Instructions - DOG

generatio


Animal Trust Center
Heidelberg - Germany
www.generatio.com

Mucous membrane swab


Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organization, follow the applicable regulations

Completed tube with mucous membrane swab

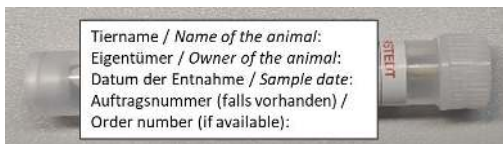

Prepare sample shipment

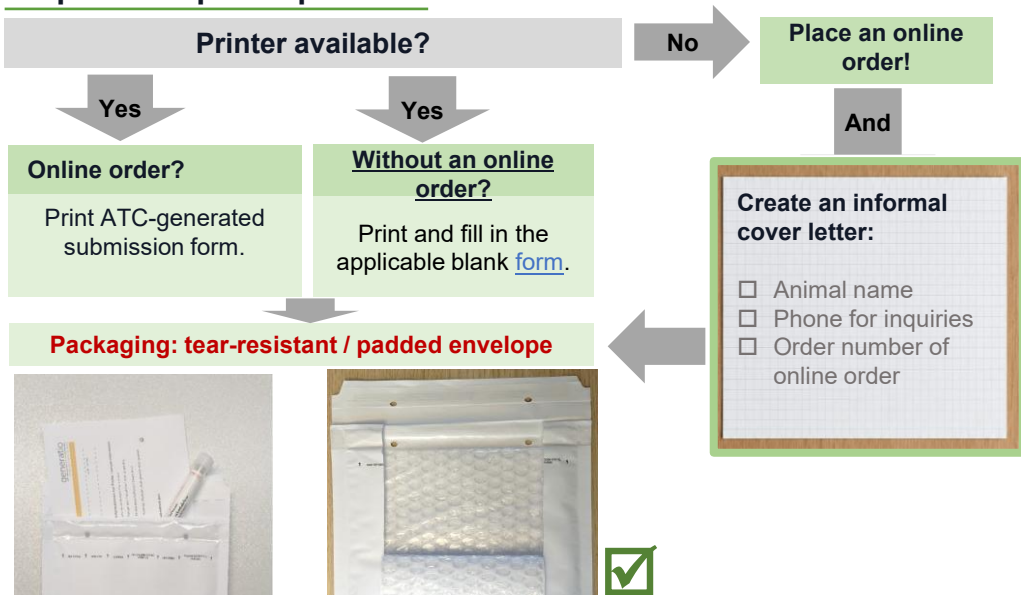

Shipping instructions



Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax

Sampling Instructions - DOG


EDTA blood sample (vet. personnel only)

Important notes:

1. *For officially recognised tests, sampling must be performed by a veterinarian.*
2. *When submitting samples for a breeding organization, follow the applicable regulations*

Preparation

☐ Label tube for EDTA blood

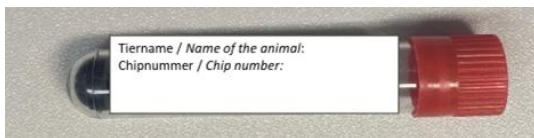

Sampling

- ➡ Draw 2–3 ml of blood.
- ➡ Carefully mix blood with EDTA by gently swirling several times.
- ➡ Please do not additionally seal the EDTA blood tube with Parafilm.

Notes

- 👉 Protect sample from direct sunlight.
- 👉 If shipping later than 2 days:
store sample tube in the refrigerator.

Sampling Instructions - DOG

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

EDTA blood sample (vet. personnel only)

Important notes:


1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organization, follow the applicable regulations

Completed EDTA blood tube



Pack the labelled EDTA blood in a break-resistant transport sleeve.

Prepare sample shipment

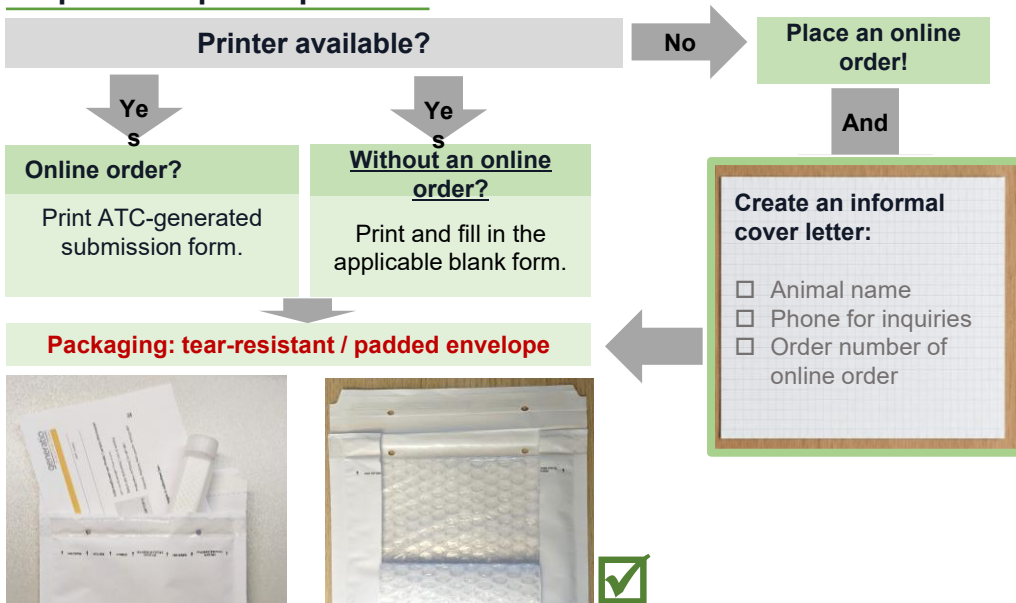

Shipping instructions



Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax